


ROSA FAST5 DON Quantitative Test for Feed and Grain

5
minutes

Product Overview

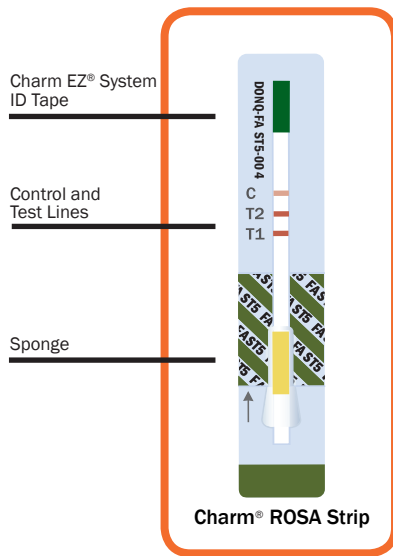

The **ROSA® FAST5 DON Quantitative Test (DONQ-FAST5)** delivers rapid and accurate results using ROSA (Rapid One Step Assay) lateral flow technology. The diluted extract is applied to the test strip, incubated for 5 minutes, and read in the Charm EZ-M system.

DONQ-FAST5 Test:

- Range: 0-1.5 ppm, 1-6 ppm, >5ppm
- Results in 5 minutes
- Negative samples may be disposed as normal waste

Run on the Charm EZ-M System:

- High volume output with batching using multiple quad incubators
- Standard curve built into reader reducing testing cost by eliminating end-user calibration
- Printable or downloadable results to any computer or LIMS system

Commodities

Charm Validated Commodities				
Barley	Brewer's Rice	Buckwheat	Corn	Corn Bran
Corn Germ Meal	Corn Gluten Feed	Corn Gluten Meal	DDGS	Hominy
Malted Barley	Milled Rice	Oats	Palm Kernal Meal	Rapeseed Meal
Rice Bran	Rough Rice	Rye	Sorghum	Soybean Meal
Triticale	Wheat	Wheat Bran	Wheat Flour	Wheat Midds
Wheat Red Dog	Wheat Screenings			



Simple Procedure



Add Sample.

Incubate.

Read results in less than 5 seconds!



For batch processing, incubate multiple strips simultaneously, then read results in the Charm EZ-M system.

Ordering Info

Order Codes	Each Kit Includes
LF-DONQ-FAST5-20K LF-DONQ-FAST5-100K LF-DONQ-FAST5-500K	Test Strips, Positive Control, Dilution Buffer

See Operator's Manual for additional kits and order codes