


# ROSA FAST5 Zearalenone Quantitative Test for Feed and Grain

5  
minutes

## Product Overview

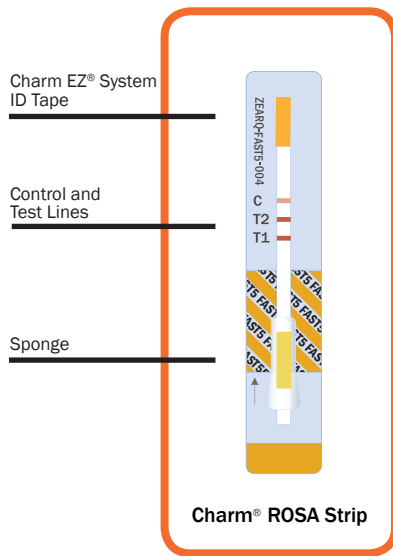

The **ROSA® FAST5 Zearalenone Quantitative Test (ZEARQ-FAST5)** delivers rapid and accurate results using ROSA (Rapid One Step Assay) lateral flow technology. The diluted extract is applied to the test strip, incubated for 5 minutes, and read in the Charm EZ-M system.

### Benefits of ZEARQ-FAST5 Test Include:

- Range: 0 to 350 ppb, 300 to 1400 ppb
- Results in 5 minutes

### Run on the Charm EZ-M System:

- High volume output with batching using multiple quad incubators
- Standard curve built into reader reducing testing cost by eliminating end-user calibration
- Printable or downloadable results to any computer or LIMS system

## Commodities

Charm® Validated Commodities				
Barley	Brewer's Rice	Corn	DDGS	Flaking Corn Grits
Milled Rice	Oats	Rice	Rough Rice	Rye
Sorghum	Soybean Meal	Wheat	Wheat Bran	Wheat Flour



## Simple Procedure



**Add Sample.**

**Incubate.**

**Read results in less than 5 seconds!**



**For batch processing, incubate multiple strips simultaneously, then read results in the Charm EZ-M system.**

## Ordering Info

Order Codes	Each Kit Includes
<b>LF-ZEARQ-FAST5-20K</b> <b>LF-ZEARQ-FAST5-100K</b> <b>LF-ZEARQ-FAST5-500K</b>	Test Strips, Positive Control, Dilution Buffer

See Operator's Manual for additional kits and order codes