



# Data Management System

## New Era of Document Control and the Power to Affect Change

The Charm® novaLINK software solves the largest data management problems facing processors and producers today. It gives you the ability to affect change directly and increase efficiency while helping lower chemical costs with the press of a button.



### Endless Possibilities

- **Scheduled Automatic Email Reports** – to be your eyes and ears via remote visibility
- **Wireless download of data and upload of test plans** – eliminating the need for connector cables (2.4 or 5.0 GHz)
- **SMART Plan Creation** – ability to create risk-based plans for a more robust quality system
- **Real-time Corporate Reporting** – results from multiple factories into one database, consolidated to a single report and can be sent out automatically.

TIMESTAMP	PLAN	LOCATION	SURFACE	RESULT RLU	PASS/FAIL	TEST NUMBER	OPERATOR	NOTES
6/23/2021 8:25:19 AM	CRUMB PREOP	MIXER PADDLES	PLASTIC	0	PASS	3	JOHN	
6/23/2021 9:19:42 AM	LINE 2	FILL VALVE 1	STAINLESS STEEL	0	PASS	4	JOHN	SUPERVISOR CONTACTED
6/25/2021 9:42:22 AM	SYRUP ROOM	MIX TANK GASKET	RUBBER	0	PASS	4	JOHN	SUPERVISOR CONTACTED
6/28/2021 10:50:58 AM	SYRUP ROOM	FILLER NOZZLE	RUBBER	44792	FAIL	3	JANE	RE-WASH & SANITIZE



The novaLINK software also adds **real-time analysis to help you see high failure areas or the need for staff retraining.**

The software can collect data in real-time from multiple facilities via wifi sync.

### Some Highlights Include:

- novaLUM II-X system includes novaLINK software with no subscription fees for local installs.
- novaLINK software only needs one (1) SQL instance and works with SQL express.
- novaLINK software economically manages multiple plants and users.
- The software can be loaded onto multiple computers with no license fees.
- Scheduled reports can automatically be sent to a designated list.
- novaLINK software hosting is available.

[www.charm.com](http://www.charm.com)

T +1.978.687.9200 | F +1.978.687.9216 | E [info@charm.com](mailto:info@charm.com) | 659 Andover Street, Lawrence, MA 01843




# How it Works

## Simple and Intuitive

The novaLINK software SMART PLAN sorts test locations into daily swabbing plans based on risk prioritization and wirelessly syncs to the novaLUM II-X system. This innovative and fully customizable software provides easy local or network installation with automatic report generation.

Customized reports with your company name and address

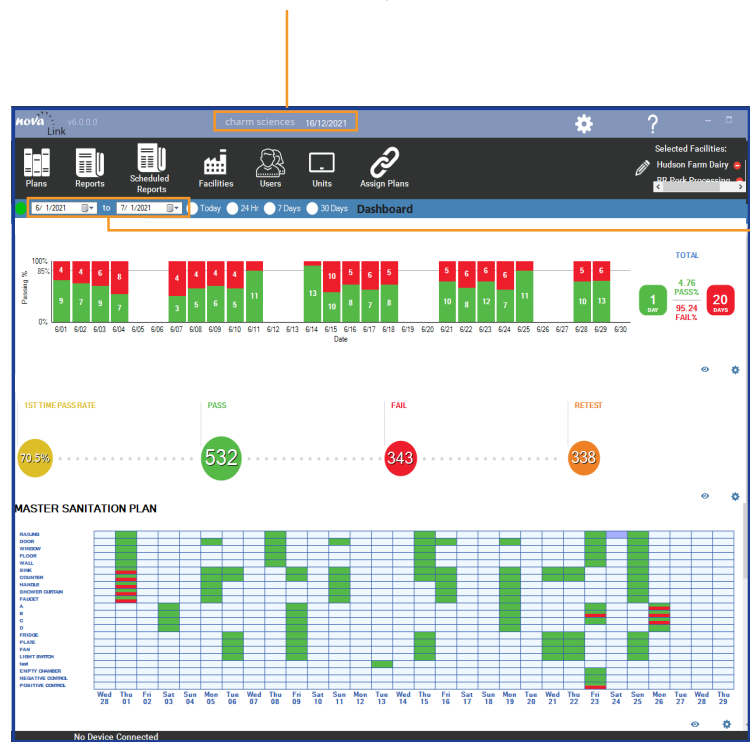


Customizable dashboard per operator

Filter by date and time

## novaLINK Dashboard

Set up multiple dashboard graphs to easily identify repetitive failures, 1<sup>st</sup>-time pass rates, and corrective actions to meet 3<sup>rd</sup> party auditing needs.

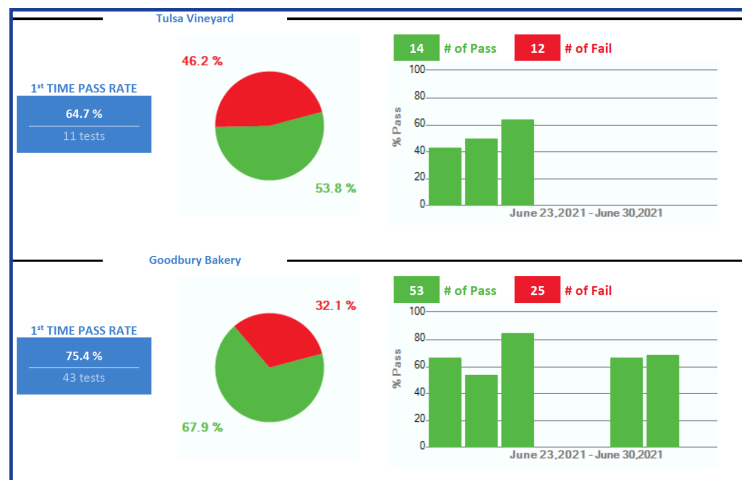

## Corporate Report

### Scheduled Email Reports

Design multiple email reports to provide key factory metrics daily to local managers and weekly to regional managers. With 6 customizable report sections novaLINK software seamlessly emails the right data to the right people.

[www.charm.com](http://www.charm.com)

T +1.978.687.9200 | F +1.978.687.9216 | E [info@charm.com](mailto:info@charm.com) | 659 Andover Street, Lawrence, MA 01843-1032, USA  
© 2022 Charm Sciences, Inc. Charm and novaLUM are registered trademarks of Charm Sciences, Inc. | MRK-1236-002 Feb-2022

